

Turny 6-Way

SV Bruksanvisning
EN User manual
DE Gebrauchsanweisung
FR Instructions d'utilisation
ES Instrucciones de uso
IT Istruzioni per l'uso

SV


EN

DE

FR

ES

IT

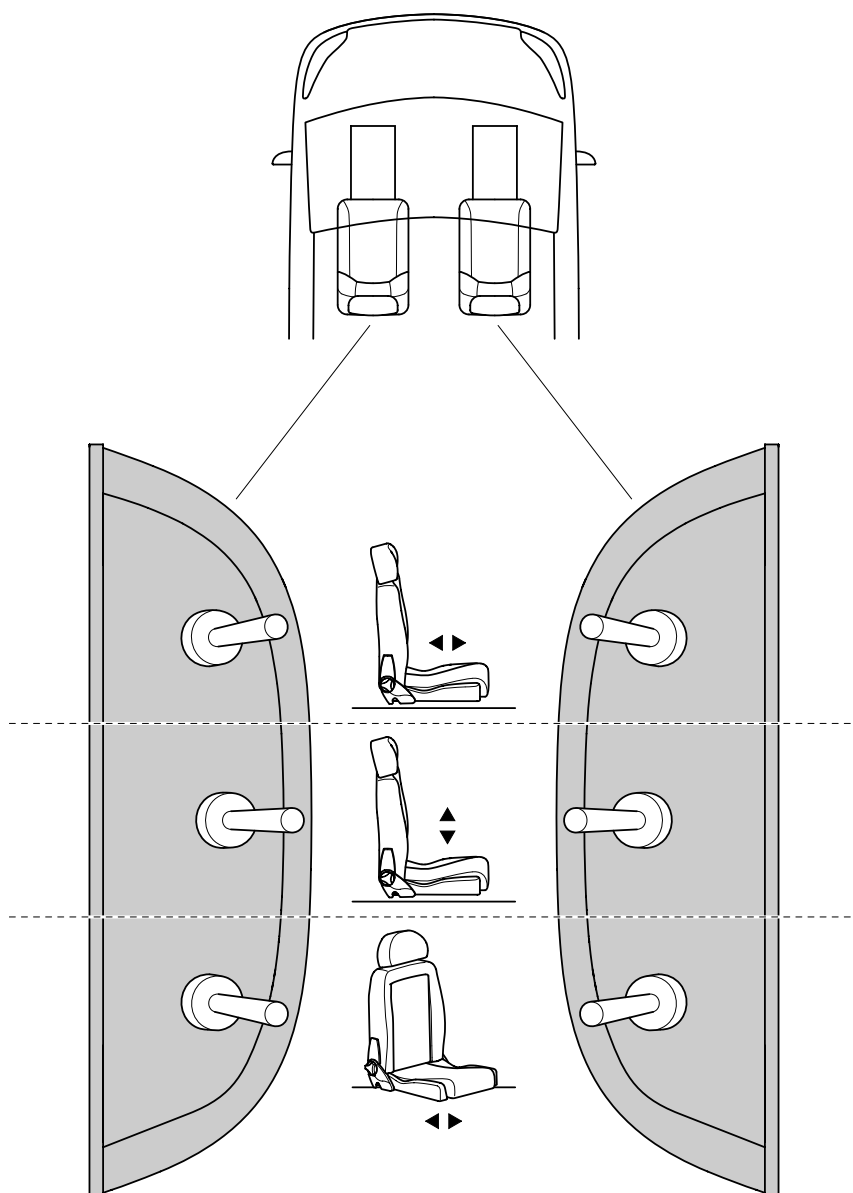






427615-2 2022/06

BraunAbility Europe AB
Åkerivägen 7
443 61 Stenkullen
Sweden

Phone: +46 302 254 00
E-mail: info@braunability.eu
www.braunability.eu

 **BraunAbility**


EN

-  *Weight capacity max 180 kg.*
-  *The Forward/Backwards and Up/Down toggle switches are always on. Ensure that they are not accidentally activated while driving. The Left/Right toggle switch is off when the vehicle's ignition is turned on.*
-  *Using the Turny 6-Way without the optional footrest causes a potential pinch risk.*
-  *Any modifications, adjustments or repairs to a Turny 6-Way must be performed by an authorised technician. Always contact your dealer for disassembly.*
-  *The control unit and the toggle switch cover shall not be used as support when transferring to and from the seat.*

Product data

Length*	1150 mm
Width	430 mm
Height	170-347 mm
Comfort adjustment	704 mm
Vertical travel	177 mm
Swivel range**	0-140° / 0-120°
Weight capacity***	180 kg
Unit weight, without seat and accessories	54 kg
Operating voltage	12V
Fuse	30A

* can be reduced 45-190 mm by installer if needed.
 ** differs between left (140°) and right (120°) side installation.
 *** refers to the total weight of the seat, its accessories and the user.

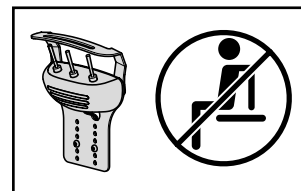

Care instructions

As necessary.....	Carefully wipe off and vacuum clean the Turny 6-Way. Do not use any liquids.
Once every two years.....	We recommend to contact your dealer for an overhaul.

Troubleshooting

Turny 6-Way is not working.....	Check that / <ul style="list-style-type: none"> - the fuse has not blown and that the vehicle's battery is charged up. - there is no obstruction blocking the selected movement. If the problem remains, contact your dealer.
Forward/backward movement is..... working slower than usually.	Probably time for re-lubrication of the aluminum profiles. Contact BraunAbility Support for further help (lubrication art no 407611).
Left/right toggle switch is working..... despite the vehicle's ignition is turned on.	Contact your dealer.
Discordant sounds are heard during operation.	Contact your dealer.